


DORMILET

PRESSURE RELIEF FOR BEDRIDDEN PEOPLE

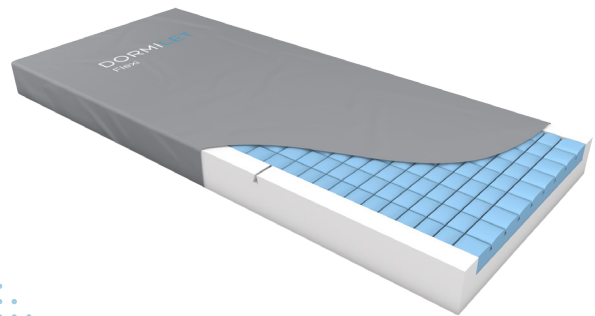

DORMILET

Flexi


MEDIUM RISK
Categories 1-2

DORMILET Flexi is the next generation of foam mattresses. It is a mattress that is designed to have the optimum effect together with the profiling functions of a bed. This combined with unique design protected cuts in the foam offers great comfort for the bedridden person combined with good pressure relieving properties.

UNIQUE STRUCTURE

DORMILET Flexi consists of two different types of foam. The bottom of the mattress and the side structure are made of firm foam. This foam frame has specially developed cutouts, which means that the mattress is aligned with the bed when for instance the back or leg support is elevated.

Inside the foam frame there is softer foam which also has a unique cutout. This cutout provides the foam with unusually great adaptability to the curvatures of the body, which offers great pressure distribution and comfort.

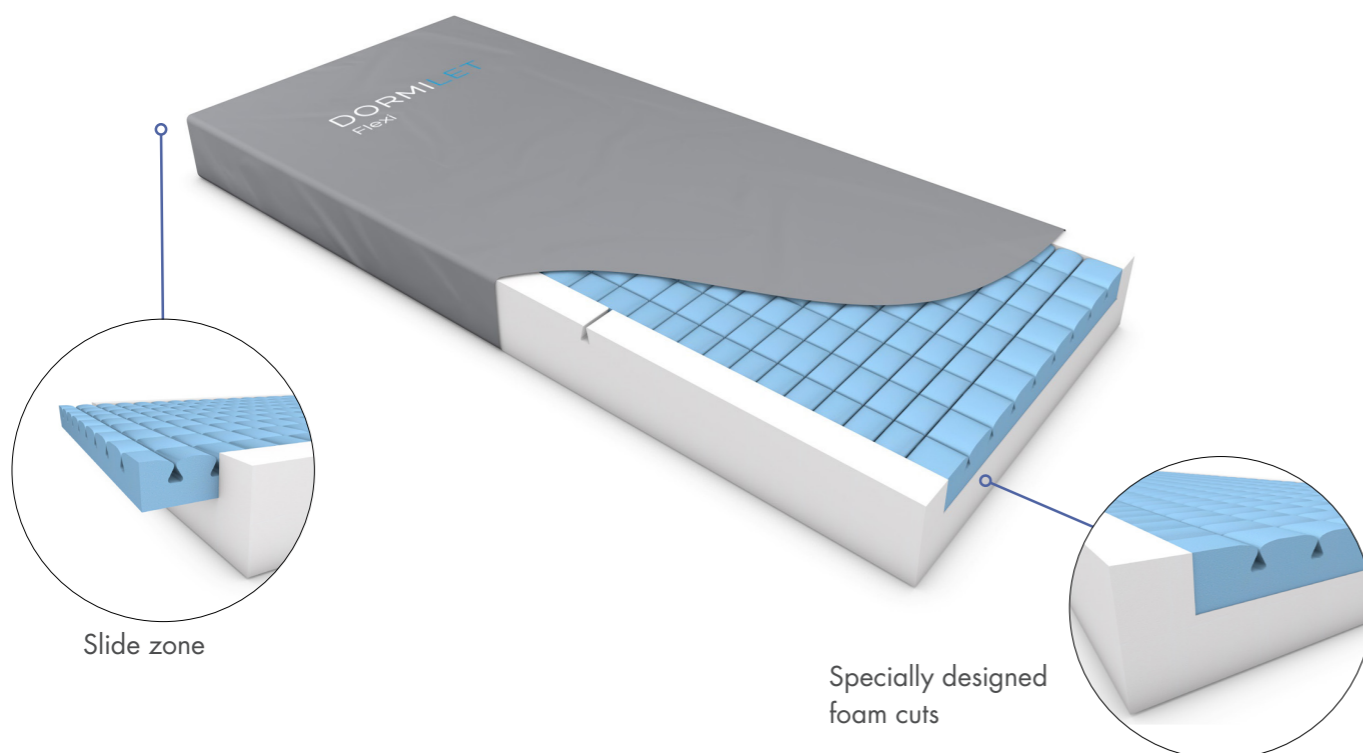

SLIDE ZONE

DORMILET Flexi has a very special slide zone. The slide zone is under the entire upper part of the body and is designed to reduce shear and friction between the bedridden person and the mattress when the back support is elevated - e.g. when a sitting position in the bed is made.

The structure means that the friction and shear normally existing between the bedridden person and the mattress is now relocated to occur inside the mattress between the two different types of foam. This not only reduces shear and friction under the back but also the pressure under the pelvis is reduced considerably.

STABILITY AND COMFORT

CMHR foam of the best quality has been used, which in combination with the special cuts and the unique slide zone gives particularly good comfort for the bedridden person.

At the same time, the stable foam frame secures comfortable and easy transfers where the user feels safe and stable - even when sitting on the edge of the mattress.

TARGET GROUP

The target group is the active/moderately active bedridden person who has up to medium risk of developing pressure sores or who has already developed sores up to and including category 2. Maximum user weight is 160 kg.

WARRANTY

There is a two year warranty on cover, foam and pump.

Item no.	Name	Size
5200000	Flexi 80	200 x 80 x 14
5200001	Flexi 85	200 x 85 x 14
5200003	Flexi 88	200 x 88 x 14

Possibility of customised mattresses in other widths. Extension pieces of 10, 15 or 20 cm can be purchased separately. See the back for additional information.



COVER AND BOTTOM

The cover of DORMILET Flexi offers multi stretch so that the foam profiles may function optimally.

Furthermore, the cover is breathable and treated with an antibacterial coating.

The seams are welded and the zipper is protected with a cover which ensures maximum protection of the foam core of the mattress. All seams are made to comply with the strictest requirements to infection control.

There are handles on the DORMILET Flexi mattresses.

CLEANING

The cover on DORMILET Flexi can be dried off with hot water with regular detergent added.

The cover can also be dried off with a 0.1 % hypochlorite solution (1,000 ppm).

In case of substantial soiling, the mattress cover first has to be cleaned with detergent and water and then a 1 % hypochlorite solution (10,000 ppm). Do not use phenols.

The cover can also be taken off the foam core and washed in the machine. The cover can be washed at up to 90° C.

The cover must air dry in a clean setting or be tumble dried at low temperature (max. 40° C).

FIRE TESTING

DORMILET Flexi has been fire tested at an independent authority.

The mattresses have undergone:

BS EN 597-1 (cigarette test)

BS EN 597-2 (match test)

BS 6807 (crib 5)

EXTENSION PIECE

For DORMILET Flexi mattresses there are extension pieces in three lengths.

10 cm
Item no.: 5100012

15 cm
Item no.: 5100013

20 cm
Item no.: 5100014